



USER MANUAL – EN

**IN 26796 Infrared LED lamp with stand inSPORTline
Elbasan**



SevenSport s.r.o. reserves the right to make any changes and improvements to its product without prior notice. Visit our website www.insportline.eu where you will find the latest version of the manual.

SAFETY INSTRUCTIONS

- Please read the manual before use and keep it for future reference.
- Do not use in places with high humidity or in water (swimming pools, saunas, etc.)
- Do not submerge in water.
- Do not use if the power cord is damaged.
- Do not cover the lamp or use it if you are tired.
- Don't look directly into the light.
- Do not use in the genital area.
- Do not use while sleeping or if you are wet.
- Do not use continuously for more than 30 min.
- Protect against impacts and falls.
- Make sure the lamp is placed on a safe and stable place.
- If you experience any discomfort during use, stop using the lamp immediately.
- If you have applied cream or oil to your skin, you need to cover that part of your body.
- Use continuously in one place for max. 30 min.
- After use, unplug. Do not pull on the power cord.
- Do not handle the product if your hands are damp or wet.
- If you clean the lamp with a damp cloth, the lamp must be dried before the next use.

WARNING: The product is not suitable for pregnant women, children, people suffering from tumors, bleeding, blood clotting disorders, burns, damaged skin or poor blood circulation.

WARNING: The product is not suitable for people sensitive to light, heat or sensitive skin.

SPECIFICATIONS

Input voltage	110-240V,50Hz/60Hz
Rated power	30 W
Light source	LED chips
Wavelength	660nm (Red light) and 850nm (Near infrared light – not visible to the human eye)

PRODUCT DESCRIPTION



- 1. LED lights
- 2. Controller
- 3. Adapter
- 4. Eye protectors

- 5. Ropes
- 6. Adjustable rod
- 7. Base

USE



- 1. Power button
- 2. Mode
- 3. Timer

- 4. Intensity
- 5. Display

1. Power switch

If the lamp is off, press the button to start the lamp. If the lamp is on, pressing the button will cause the fluorescent lamp to flash. Press it again to turn off the lamp.

2. Mode

The lamp has 3 lighting modes:

- a. 660 nm + 850 nm LED radiation
- b. 660 nm radiation
- c. 850 nm radiation

After starting, radiation is on and press the button to change the mode.

3. Timer

You can set the timer in 6 levels: 5 min, 10 min, 15 min, 20 min, 25 min, 30 min.


The default is no timer. Press the button to set the timer in a cycle.

4. Intensity

There are 5 intensity levels to choose from: L0 – L4. The setting changes in a cycle.

CABLE ATTACHMENT

Hook the carabiners to the eyelets on the lamp as shown in the picture.

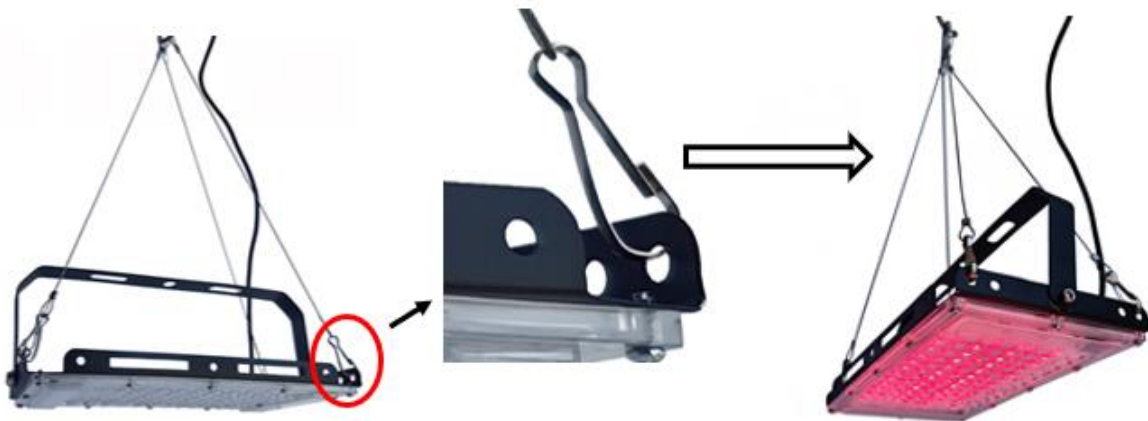

PLACING ON THE TABLE

You can adjust the angle using the support on the back.



STAND

Step 1

Remove the packing material and spread the legs of the stand.

Turn it clockwise to loosen the stand, turn it counter-clockwise to tighten it.


Step 2

Attach the lamp to the stand. Turn clockwise to tighten the lamp, turn counterclockwise to loosen.

Plug the lamp into the adapter.

Step 3

Plug the adapter into the electrical outlet. current. Choose the mode, intensity and duration of use.

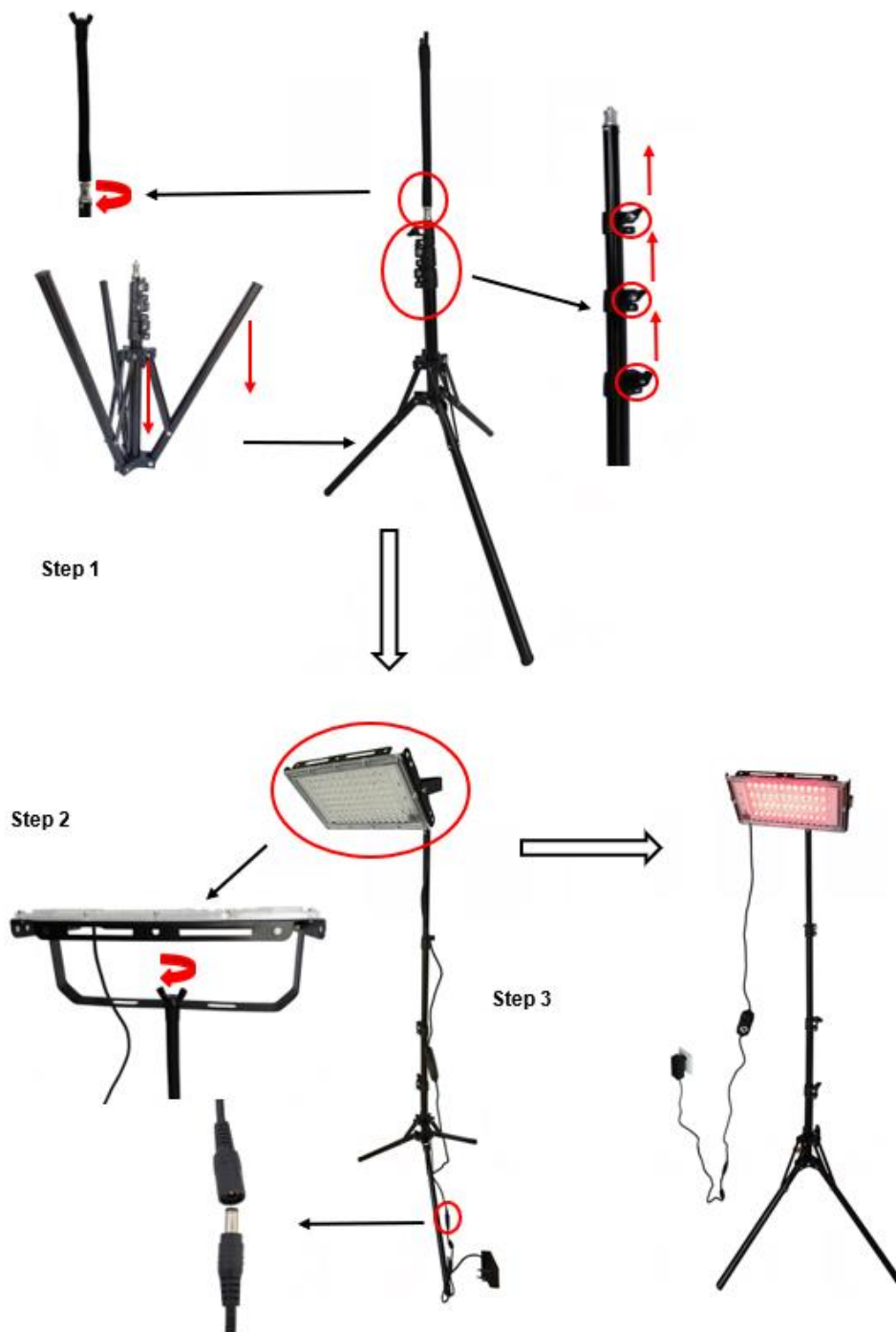

ENVIRONMENT PROTECTION

After the product lifespan expired or if the possible repairing is uneconomic, dispose it according to the local laws and environmentally friendly in the nearest scrapyard.

By proper disposal you will protect the environment and natural sources. Moreover, you can help protect human health. If you are not sure in correct disposing, ask local authorities to avoid law violation or sanctions.

Don't put the batteries among house waste but hand them in to the recycling place.

TERMS AND CONDITIONS OF WARRANTY, WARRANTY CLAIMS

General Conditions of Warranty and Definition of Terms

All Warranty Conditions stated here under determine Warranty Coverage and Warranty Claim Procedure. Conditions of Warranty and Warranty Claims are governed by Act No. 89/2012 Coll. Civil Code, and Act No. 634/1992 Coll., Consumer Protection, as amended, also in cases that are not specified by these Warranty rules.

The seller is SEVEN SPORT s.r.o. with its registered office in Strakonická street 1151/2c, Prague 150 00, Company Registration Number: 26847264, registered in the Trade Register at Regional Court in Prague, Section C, Insert No. 116888.

According to valid legal regulations it depends whether the Buyer is the End Customer or not.

“The Buyer who is the End Customer” or simply the “End Customer” is the legal entity that does not conclude and execute the Contract in order to run or promote his own trade or business activities.

“The Buyer who is not the End Customer” is a Businessman that buys Goods or uses services for the purpose of using the Goods or services for his own business activities. The Buyer conforms to the General Purchase Agreement and business conditions.

These Conditions of Warranty and Warranty Claims are an integral part of every Purchase Agreement made between the Seller and the Buyer. All Warranty Conditions are valid and binding, unless otherwise specified in the Purchase Agreement, in the Amendment to this Contract or in another written agreement.

Warranty Conditions

Warranty Period

The Seller provides the Buyer a 24 months Warranty for Goods Quality, unless otherwise specified in the Certificate of Warranty, Invoice, Bill of Delivery or other documents related to the Goods. The legal warranty period provided to the Consumer is not affected.

By the Warranty for Goods Quality, the Seller guarantees that the delivered Goods shall be, for a certain period of time, suitable for regular or contracted use, and that the Goods shall maintain its regular or contracted features.

Batteries

6-month battery warranty – we guarantee that battery's nominal capacity does not fall below 70% of its total capacity within 6 months of the product's sale.

The Warranty does not cover defects resulting from (if applicable):

- User's fault, i.e. product damage caused by unqualified repair work, improper assembly, insufficient insertion of seat post into frame, insufficient tightening of pedals and cranks
- Improper maintenance
- Mechanical damages
- Regular use (e.g. wearing out of rubber and plastic parts, moving mechanisms, joints, wear of brake pads/blocks, chain, tires, cassette/multi wheel etc.)
- Unavoidable event, natural disaster
- Adjustments made by unqualified person

- Improper maintenance, improper placement, damages caused by low or high temperature, water, inappropriate pressure, shocks, intentional changes in design or construction etc.

Warranty Claim Procedure

The Buyer is obliged to check the Goods delivered by the Seller immediately after taking the responsibility for the Goods and its damages, i.e. immediately after its delivery. The Buyer must check the Goods so that he discovers all the defects that can be discovered by such check.

When making a Warranty Claim the Buyer is obliged, on request of the Seller, to prove the purchase and validity of the claim by the Invoice or Bill of Delivery that includes the product's serial number, or eventually by the documents without the serial number. If the Buyer does not prove the validity of the Warranty Claim by these documents, the Seller has the right to reject the Warranty Claim.

If the Buyer gives notice of a defect that is not covered by the Warranty (e.g. in the case that the Warranty Conditions were not fulfilled or in the case of reporting the defect by mistake etc.), the Seller is eligible to require a compensation for all the costs arising from the repair. The cost shall be calculated according to the valid price list of services and transport costs.

If the Seller finds out (by testing) that the product is not damaged, the Warranty Claim is not accepted. The Seller reserves the right to claim a compensation for costs arising from the false Warranty Claim.

In case the Buyer makes a claim about the Goods that is legally covered by the Warranty provided by the Seller, the Seller shall fix the reported defects by means of repair or by the exchange of the damaged part or product for a new one. Based on the agreement of the Buyer, the Seller has the right to exchange the defected Goods for a fully compatible Goods of the same or better technical characteristics. The Seller is entitled to choose the form of the Warranty Claim Procedures described in this paragraph.

The Seller shall settle the Warranty Claim within 30 days after the delivery of the defective Goods, unless a longer period has been agreed upon. The day when the repaired or exchanged Goods is handed over to the Buyer is considered to be the day of the Warranty Claim settlement. When the Seller is not able to settle the Warranty Claim within the agreed period due to the specific nature of the Goods defect, he and the Buyer shall make an agreement about an alternative solution. In case such agreement is not made, the Seller is obliged to provide the Buyer with a financial compensation in the form of a refund.

CZ SEVEN SPORT s.r.o.

Registered Office: Strakonická 1151/2c, Praha 5, 150 00, ČR
Headquarters: Dělnická 957, Vítkov, 749 01
Warranty & Service: Čermenská 486, Vítkov 749 01

CRN: 26847264
VAT ID: CZ26847264
Phone: +420 556 300 970
E-mail: eshop@insportline.cz
reklamace@insportline.cz
servis@insportline.cz

Web: www.inSPORTline.cz

SK Stores inSPORTline SK, s.r.o.

Headquarters, warranty & service center: Električná 6471, Trenčín 911 01, SK

CRN: 46259317
VAT ID: SK2023299729
Phone: +421(0) 917 700 098
E-mail: objednavky@insportline.sk
reklamacie@insportline.sk
servis@insportline.sk

Web: www.inSPORTline.sk

About shipping

